

## Importing Airtraq video files into eMeds Elite.

The new Airtraq A-390 Wi-Fi Camera allows EMS clinicians to upload videos of intubations into eMeds Elite. Clinicians are to upload ALL successful AND attempted intubations performed with the Airtraq into eMEDS Elite. This document only covers technical aspects of the Airtraq A-390



### Airtraq A-390 Connections and Charging

A USB-A to USB-Micro Cable will be supplied with each Airtraq A-390 so that clinicians can download videos and charge the camera. The Airtraq can be charged on apparatus using the charger for the Sonim apparatus cell phones. The Airtraq A-390 can be turned on to check the battery by depressing the POWER ON switch (pictured below). The Airtraq A-390 can operate for 240 minutes with the Wi-Fi Disabled and requires 110 minutes to charge.




### Recording Videos and Images

When the Airtraq is turned on it should start recording immediately. The screen will display STOP REC and should record for 10 minutes. If the STOP REC is not displayed press the REC button prior to the attempt. To capture a still image, simply tap the screen to save a picture.

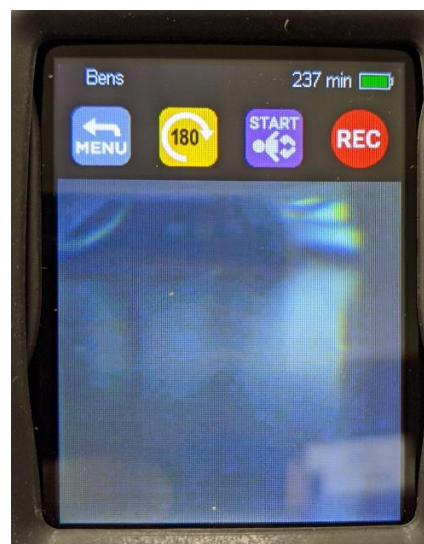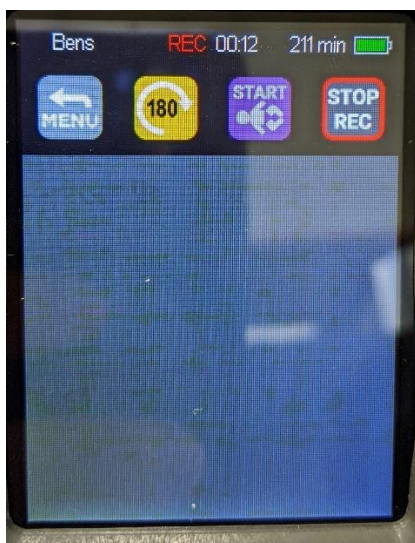


### Menu Screen

MCFRS does not use the Wi-Fi in operations and turning the Wi-Fi on disables device recording, so **leave the Wi-Fi disabled**.

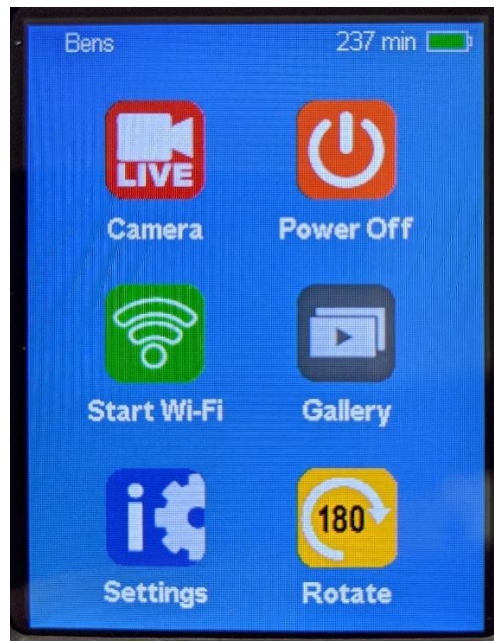


In settings, Auto Recording should default to ON, and REC Time Limit should be 10 min.

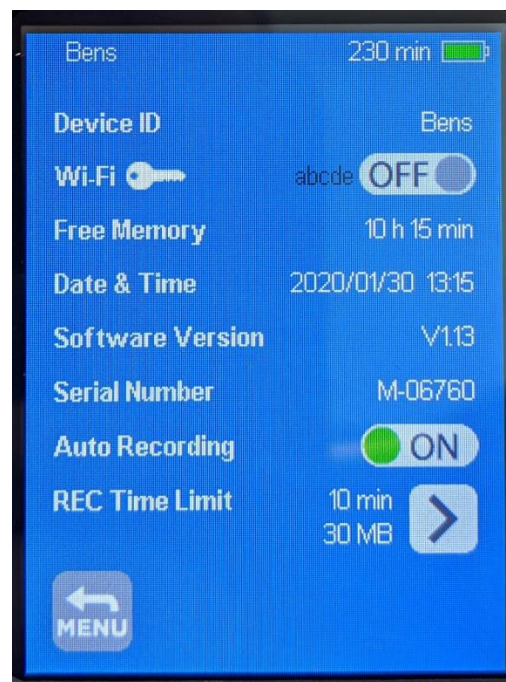

### Upload Files into eMeds Elite

Connect The Airtraq A-390 to your ePCR with the supplied USB-A to USB-Micro cable.

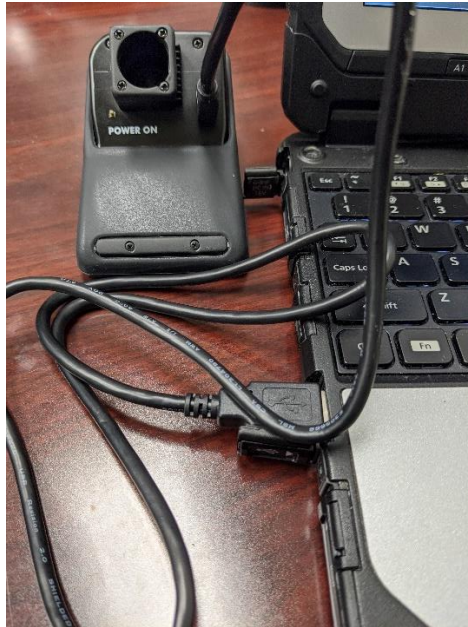


In your ePCR report:

1. Click Menu




2. Click Attachment

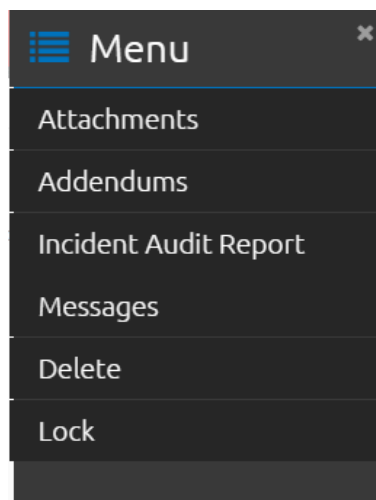


3. Click Add New




4. Click Browse

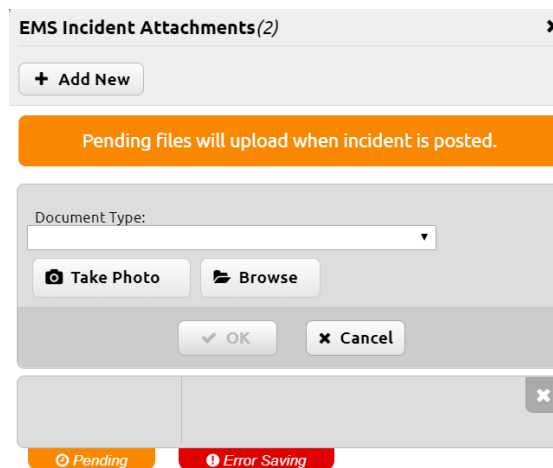

5. In the Explorer Window, under "Computer" on the Left side click "AIRTRAQ CAM". Your recorded intubation attempts will be in the "videos" folder. Double click the "videos" folder to open it.

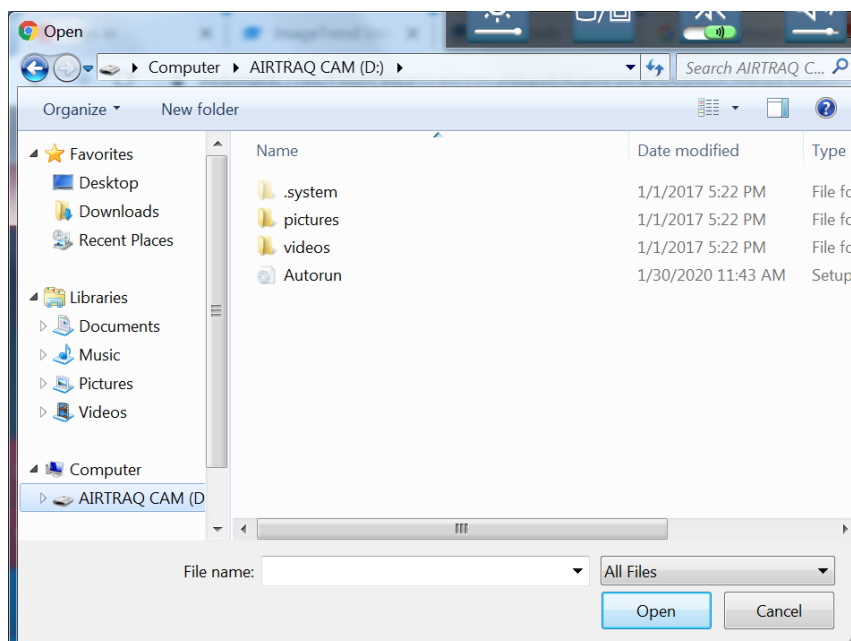

- The filenames have the date and time of every video created. Click the oldest file created during your incident, then click open, then click the upper OK

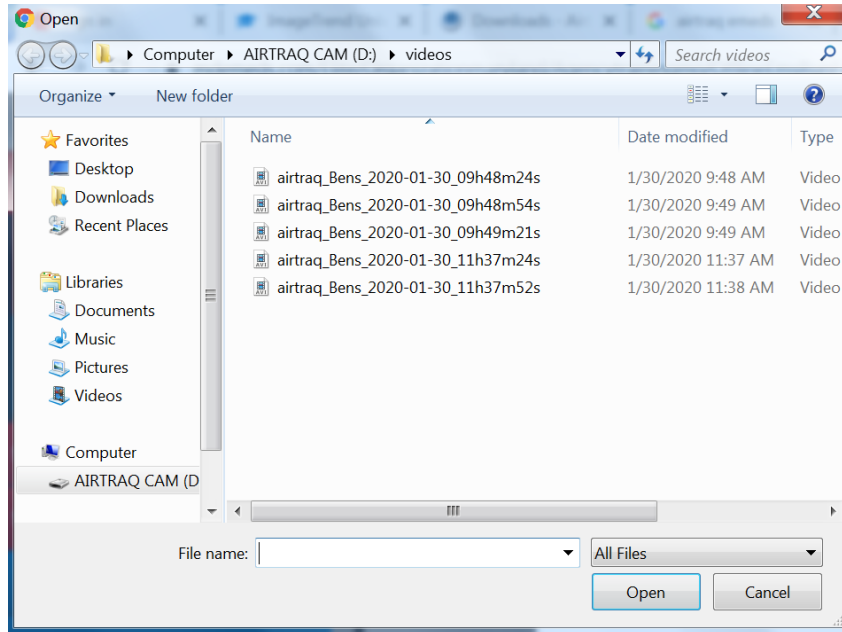

- Repeat steps 3-6 for every video created during your incident.
- Click OK



9. Complete your report as normal. Be sure to add a procedure for “Airway – Laryngoscopy, Video”. Note: It will take a longer time to post the incident.

### Procedure

+ Add Another

✓ OK

✕ Cancel

Prior to Arrival?:

No

Yes

⊖

Role/Type of Person Performing the Procedure:

✕ Paramedic

▼

☰

⊖

Date/Time:

01/30/2020

13:32:33

🕒

Performed By:

✕ PRACTICE, ALS (15000018)

▼

☰

⊖

Procedure Name:

✕ Airway - Laryngoscopy, Video

▼

☰

⊖